

PART OF THE **DURAN®** BOTTLE SYSTEM

EFFICIENT

DURAN® HPLC RESERVOIR BOTTLE



DWK LIFE SCIENCES
DURAN WHEATON KIMBLE

DURAN® HPLC RESERVOIR BOTTLE CLEAR



conical base, with GL 45



The DURAN® reservoir bottles have been specially designed for the storage and delivery of high-performance liquid chromatography (HPLC) mobile phases and solvents. Reservoirs have a one central cavity for the inlet filter to allow the most efficient delivery of mobile phase to HPLC analyzers. A safer and more convenient alternative to flat bottomed glass media bottles, the conical bottom reservoirs allow a greater recovery of solvent, without having to tilt reservoir and risk a solvent spill. The larger reservoirs allow longer, uninterrupted analyses. These are especially useful during a weekend or lengthy chromatographic analysis when there might be a risk of running out of solvents during unattended operation.


Cat. No.	Capacity (mL)	DIN Thread (GL)	d (OD) (mm)	h (mm)	Pack Unit
21 841 54 07	1 000	45	110	275	1
21 841 63 09	2 000	45	145	313	1
21 841 69 09	3 500	45	165	353	1
21 841 73 05	5 000	45	200	400	1
21 841 86 01	10 000	45	250	482	1

Cat. No.	Description	Pack Unit
Accessories		
29 231 28 01	DURAN® 3-port HPLC GL 45 delivery cap 1/4" - 28 thread for 3.2 mm OD tubing PTFE Complete with three 1/4" - 28 ETFE plugs and a TFE / propylene O-ring	1
29 231 00 03	DURAN® ETFE sealing plug 1/4" - 28 thread for 3-port HPLC GL 45 delivery cap	1

DURAN® GL 45 HYDRA SQUARE BOTTLE WITH FOUR GL 45 SIDE NECKS



HPLC purging bottle 1000 ml with five GL 45 inlets

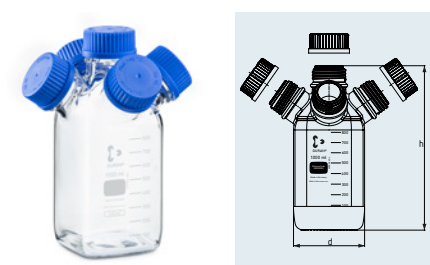

DURAN® bottles are ideal for the storage and delivery of high-performance liquid chromatography (HPLC) mobile phases and solvents. Increasing the number of GL 45 inlets on the bottles adds convenience: The five GL 45 necks on a purging bottle increases safety and reduces the risk of cross contamination and air intake.


Cat. No.	Capacity (mL)	DIN Thread (GL)	d (OD) (mm)	h (mm)	Pack Unit
21 850 54 53	1 000	45	94	222.5	1



DURAN® GL 45 BOTTLE 2000 ML WITH GL 45 SIDE NECK



HPLC Reservoir Bottle 2000 ml with two GL 45 inlets

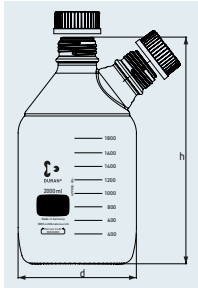

DURAN® bottles are ideal for the storage and delivery of high-performance liquid chromatography (HPLC) mobile phases and solvents. Increasing the number of GL 45 inlets on the bottles adds convenience: One additional side neck allows additional solvent to be added without disturbing the main delivery cap.

Cat. No.	Capacity (mL)	DIN Thread (GL)	d (OD) (mm)	h (mm)	Pack Unit
21 990 07 01	2 000	45	136	265.5	1

